

AGK 2000



Solid particle aerosols produced from suspensions, e.g., with biological agents, and solutions such as NaCl and KCl, are required for numerous research, development, and quality assurance applications. NaCl/KCl aerosols or aerosols produced from biological agents are prescribed as test aerosols in various standards to ensure the comparability of filter media, measuring equipment, and filters. To be called such, test aerosols must be generated consistently regarding particle size distribution and particle concentration over the test period. Furthermore, it must be possible to reproducibly produce the particle size distribution and the concentration. A specially developed nozzle ensures that these requirements are met by preventing the crystallization of the salt crystals at the nozzle outlet. The particle size spectrum can thus be adjusted reproducibly in the range of approx. 5 nm up to 15 μm , depending on the concentration of the solution.

OPERATION PRINCIPLE

GENERATION OF SOLID PARTICLES OUT OF SUSPENSIONS, SOLUTIONS, AND BIOLOGICAL AGENTS

Two-substance nozzles are mainly used for dispersing liquids, suspensions, and solutions. For dispersing salt solutions, the conventional atomization methods, such as that of the Collison Nebulizer, which was also built by Palas®, are not particularly suitable since salt crystals precipitate at the nozzle outlet and lead to partial obstruction of the nozzle system.

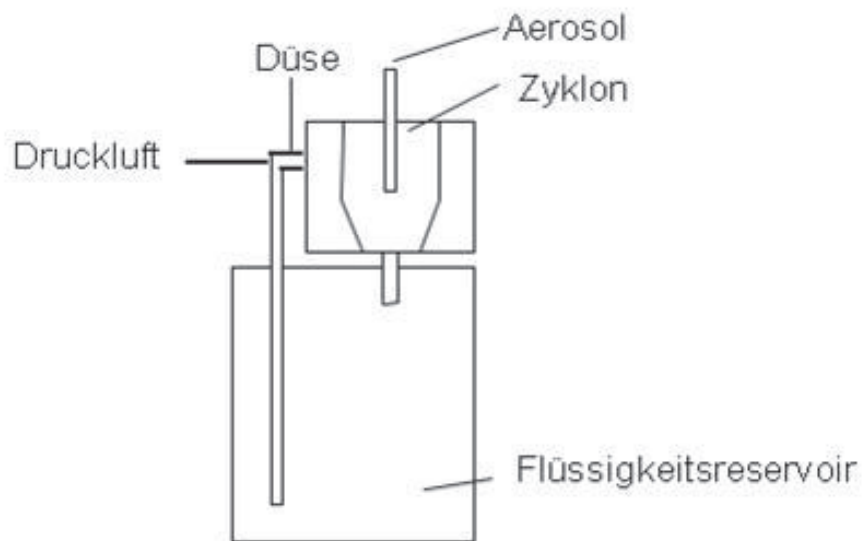


Fig. 1: Schematic diagram of the AGK 2000 with cyclone

A special nozzle developed by Palas® enables the atomization of salt solutions with the highest dosing constancy. Depending on the solution concentration of salts in the reservoir, the particle size of the aerosol can be influenced. Higher solution concentrations lead to larger particles.

Dependency of the particle distribution from the concentration of the solution (Figs. 2 and 3). Measuring device: welas® digital system from Palas®.

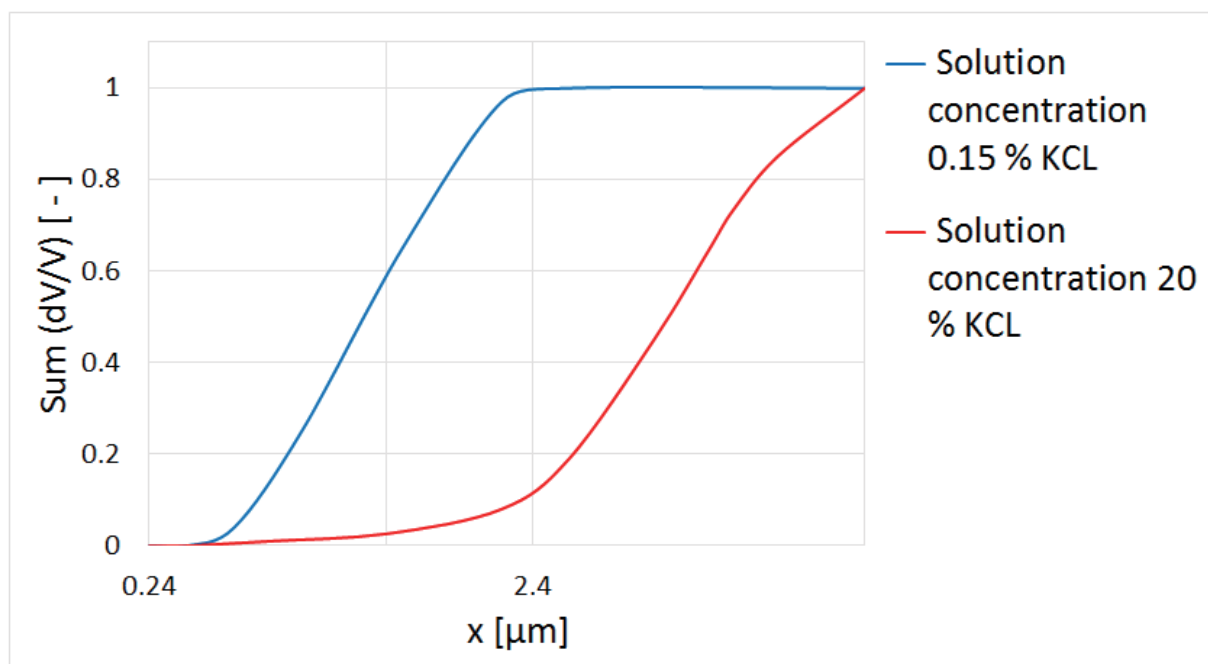


Fig. 2: Representation of the volume-based integral distribution of a 20 % KCl solution and a 0.15 KCl solution

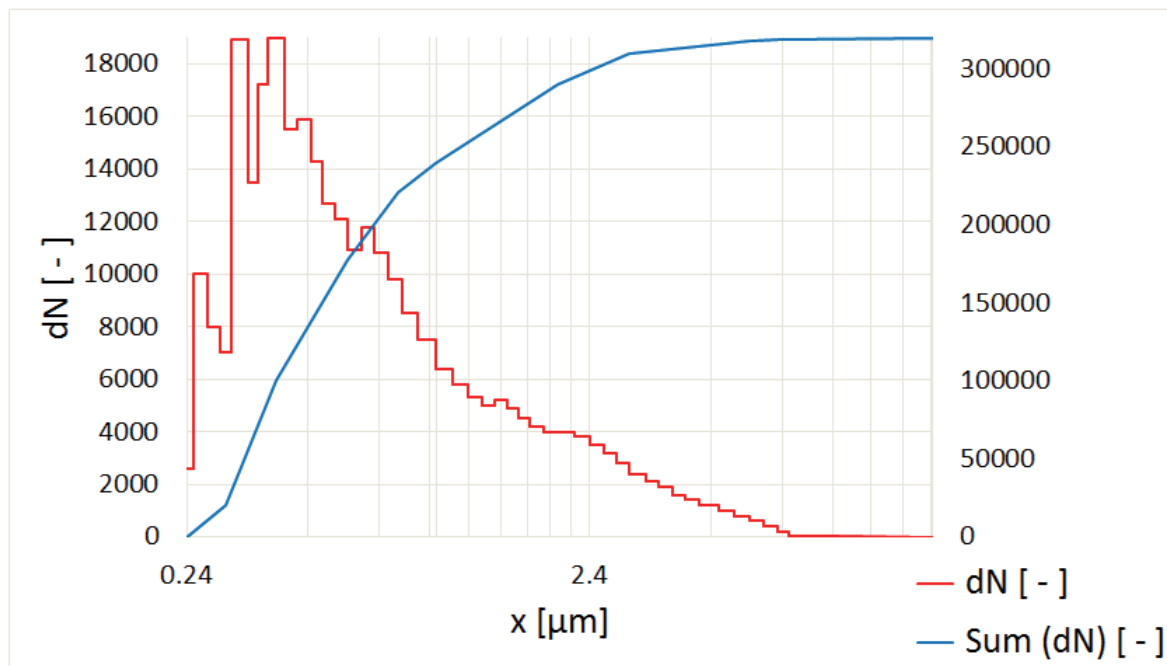


Fig. 3: Representation of the number size distribution of a 20 % KCl solution

Extensions/Accessories

The AGK 2000 has a straight or curved drying column and is optionally pressure-resistant up to 10 bar.

BENEFITS

- Excellent short-term and long-term dosing constancy
- Wide adjustable particle size range
- Easy filling of the reservoir
- Large reservoir (500 cm³)
- Robust design, proven in industrial applications
- Easy to operate
- Reliable function, high reproducibility
- Little maintenance required
- Reduces your operating costs

DATASHEET

Volume flow	3 – 10 l/min
Weight	Approx. 3 kg
Particle material	NaCl, KCL, biological agents and other particles in suspensions
Dosing time	Several hours nonstop
Maximum particle number concentration	Ca. 10^7 particles/cm ³
Particle size range	0.005 – 15 µm
Carrier/dispersion gas	Random (generally air)
Pre-pressure	4 – 8 bar
Compressed air connection	Quick coupling
Aerosol outlet connection	Ø _{inside} = 20 mm, Ø _{outside} = 30 mm
Volume flow (accessories)	6 – 20 l/min (drying column)
Filling quantity	300 ml

Applications

- Filter industry:
 - Car interior filters
 - ASHRAE room air filters
 - Engine air filters
 - Respiratory filters
- Chemical and pharmaceutical industry
- Generation of tracer particles
- Flow visualization
- Aerosol research

As an aerosol technology expert, Palas® Germany is committed to providing users with solutions for the generation, conditioning, measurement and analysis of aerosol particles. Based on the unique advantages of its own technology, Palas® developed a variety of application cases in ambient air quality monitoring, particle filtration performance testing and various scientific research fields. Palas Instruments (Shanghai) Co., Ltd. is a wholly owned subsidiary of Hong Kong Palas (Asia) Limited. As one of the global branches of Palas GmbH, it has legally obtained the Palas trademark authorized by Palas GmbH in Exclusive use rights in China and Asia.

As a company that has passed the ISO 9001:2015 quality management system certification, Palas®'s test rig solutions can execute particle filtration performance tests for filter media and filter elements according to applicable international, national and regional standards. In terms of environmental protection, Palas®'s equipment meets the requirements of multiple environmental monitoring standards (EN 15267, EN 16450, HJ653, GBZ/T 192.6, etc.) for indoor and ambient PM2.5, PM10, particle number size distribution monitoring and analysis.

Palas (Asia) Limited, Hong Kong

Operational Office Address:

Palas Instruments (Shanghai) Co., Ltd.

5th Floor, Building 6C, No. 650 Shunqing Rd, Song Jiang District, 201612 Shanghai

Hotline: +86 400 784 6669

Email: info@palas.com.hk

Website: www.palas.de/en