



RBG SYSTEM

AEROSOL GENERATOR

For dust and powders

Made in Germany

Fast. Simple. Reproducible.

RBG SYSTEM

Low-concentration solid particle aerosols produced from dust and powders are required for many applications in research, development, and quality assurance and for the calibration of particle measurement devices.

A high level of dispersion constancy is required for those substances to the point of lowest dosing rates. Additionally, a good reproducibility during aerosol generation must be guaranteed.



The **RBG SYSTEM** fulfills these requirements for mass flows between 0.04 up to 800 g/h.

A special advantage:

The dosing and dispersion systems are fast and simple adaptable to different application conditions.

This allows for example:

- Use of different powders - dusts
- Use of different carrier gases
- Overpressure operation up to 10 bar
- Negative pressure operation down to 600 mbar absolute

Application examples



FILTER TESTING




**DEVELOPMENT OF NEW FILTER
AND FILTER MATERIAL**



CALIBRATION



PARTICLE IMAGE VELOCIMETRY

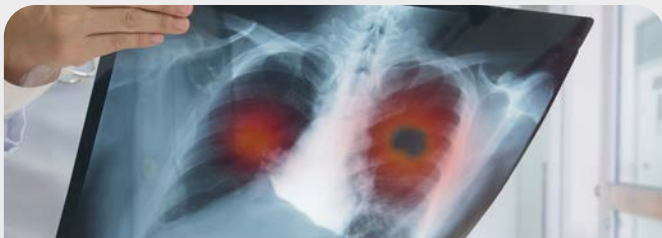

INHALATION TEST



MASK CHECK


Principle of operation

The **RBG SYSTEM** from Palas GmbH is characterized by its precision and easy handling.

The material inside the reservoir is moved upwards at a defined speed by a stepper motor and is collected by the rotating brush which transports the powder into the dispersion head. A high velocity gas stream blows the powder out of the brush and disperses it in the dispersion nozzle. The mass flow can thus be easily and reproducibly determined from the set feed rate and the cross-section of the transport piston, as well as the packing density of the powder.

The dispersion gas flow can be set manually or fully automated with the integrated air flow control. Four different dispersing covers can be used for optimal dispersion.

The **RBG SYSTEM** can be set manually via the display keys or via a connected computer with the Palas® control software included in the shipment.

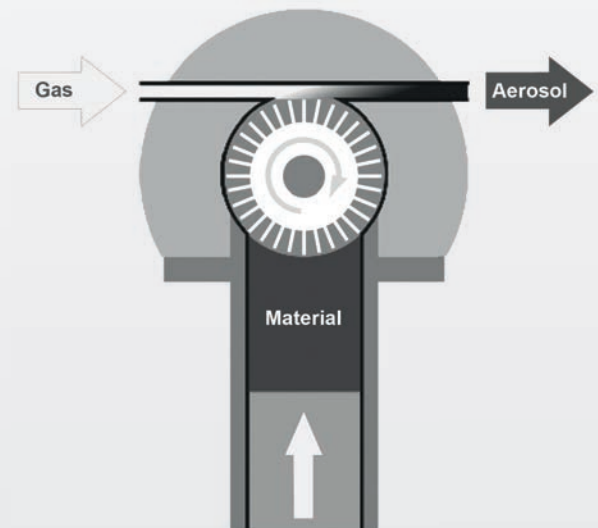

RBG SYSTEM

	RBG PROFESSIONAL	RBG BASIC	RBG SOLO
Feed stock reservoir with up to 110 mm filling height	X	X	X
USB connection for digital remote control	X	X	X
Operation with different carrier gases, e. g., nitrogen	X	X	X
Fast change of the reservoir in approx. 30 s	X	X	X
LCD functional display with keyboard	X	X	X
Robust plastic housing	X	X	X
Electrical supply by wide-range power supply unit (24 V)	X	X	X
Automatic volume air flow control with integrated mass flow meter & valve	X		X
Capsuled maintenance-free motor for dispersion brush, pressure-tight up to 10 bar	X		
Integrated pump for operation without compressed air supply			X

Special advantages and benefits

RELIABILITY


- High reproducibility for measurements due to very good long-term dosing consistency
- Ideal for dispersion of many non-cohesive dusts

FLEXIBILITY

- Fast adaptation to changing test conditions with different solids containers and dispersion lids
- Wide application range up to 10 bar over pressure and 0.4 bar negative pressure and with different carrier gases
- Fast response

USER FRIENDLINESS

- Easy to operate due to automatic monitoring of the volume flow and LCD display for all parameters
- Convenient remote control of the **RBG SYSTEM** via supplied Palas® software via USB interface

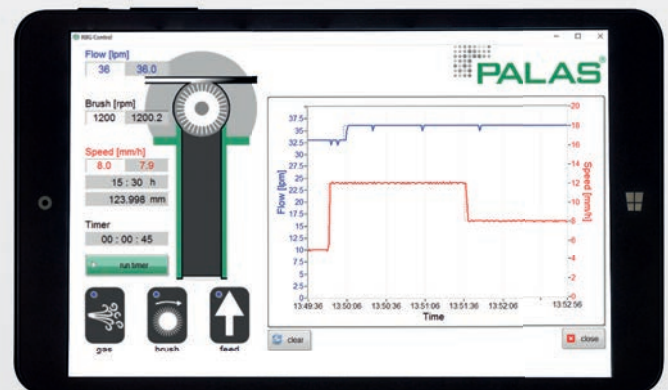

Technical features

Particle material	Non-cohesive powders and bulks
Volume flow	Depending on the dispersion cover: 8 – 180 NI/min (RBG PROFESSIONAL) 8 – 85 NI/min (RBG BASIC) 8 – 40 NI/min (RBG SOLO)
Max. particle number concentration	Approx. 10^7 particles/cm ³
Mass flow (particles)	0.04 – 800 g/h (with an assumed compacted density of 1 g/cm ³)
Particle size range	0.1 – 200 µm
Pre-pressure	4 – 13 bar (RBG PROFESSIONAL) 4 – 8 bar (RBG BASIC)
Negative pressure	0.4 bar
Maximum over pressure	10 barg (RBG PROFESSIONAL) 0.2 barg (RBG BASIC) 0.1 barg (RBG SOLO)
Filling quantity	2.7 – 88 g (with an assumed compacted density of 1 g/cm ³)
Interfaces	USB type B
Dimensions (H • W • D)	515 • 330 • 240 mm
Weight	Approx. 15 kg (RBG PROFESSIONAL, RBG BASIC) Approx. 19 kg (RBG SOLO)



Palas®, founded in 1983 and headquartered in Karlsruhe, Baden-Württemberg, Germany, specializes in the development of high-precision instruments for the generation, measurement and analysis of aerosol particles and is the world's leading developer and manufacturer in this field. Palas Instruments (Shanghai) Co., Ltd. is a wholly owned subsidiary of Hong Kong Palas (Asia) Limited. As one of the global branches of Palas GmbH, it has legally obtained the Palas trademark authorized by Palas GmbH in Exclusive use rights in China and Asia.

As an ISO 9001:2015 certified company, Palas® helps customer to achieve the testing capability of filter media and filter element in accordance with ISO 16890, ISO 11155-1, ISO 5011, ISO 29463-3, EN 1822-3, ISO 17536, ISO 12500 via tester solution. In terms of environmental protection, Palas® also complies with a number of environmental monitoring industry standards (EN 15267, HJ653, etc.).

Palas (Asia) Limited, Hong Kong

Operational Office Address:

Palas Instruments (Shanghai) Co., Ltd.

5th Floor, Building 6C, No. 650 Shunqing Rd, Song Jiang District, 201612 Shanghai

Hotline: +86 400 784 6669

Email: info@palas.com.hk

Website: www.palas.com.hk